



LABMASTER
LUCIA™

Veterinary diagnostics

**QUICK AND ACCURATE RESULTS -
WITH EASE**

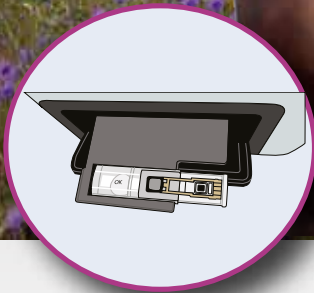

Insert the cassette
into the analyzer

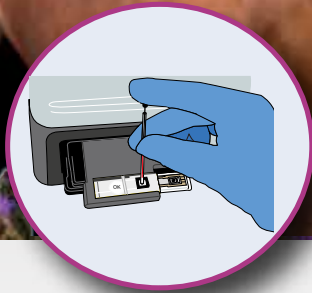

Dispense the prepared sample
into the cassette



Result can be transferred to
Windows or Android App

LABMASTER LUCIA™ VET ANALYZER

- measures analytes in whole blood
- a quantitative result in a few minutes
- a versatile platform for multiple test kits
- maintenance-free
- easy to transport

POINT-OF-CARE TEST KITS

- TEST CASSETTE with reagents and washing buffer included
- NFC CARD with test-specific parameters
- DILUTION TUBE for preparing the sample

EQUINE SAA (eSAA) TEST KIT

- for measuring the concentration of Equine Serum Amyloid A (eSAA) in whole blood

PIPELINE TESTS

- Foal IgG
- Canine CRP

FUTURE DEVELOPMENT

- Feline SAA
- Progesterone
- Testosterone



Legal disclaimer

Products designed, developed and manufactured by Labmaster Ltd. meet their specifications when stored and used according to the instructions for use. Labmaster LUCIA™ is a trademark designed and owned by Labmaster Ltd.

Products presented may not be available in all countries. Please check availability from your local representative.

LABMASTER LUCIA™

EQUINE SAA

INTENDED USE

The Labmaster LUCIA™ Equine SAA test is a veterinary diagnostic test for the quantitative determination of Equine Serum Amyloid A (eSAA) from whole blood to assess infection and inflammation status in horses.

SPECIFICATION

SAMPLE TYPE	Whole blood
SAMPLE VOLUME	10 µL
MEASUREMENT TIME	6 min
MEASURING RANGE	10–1500 mg/L
STORAGE	2–8 °C

Take advantage of the Labmaster LUCIA™ Vet Analyzer leasing programme!

We offer a full ownership of the analyzer after a 24-month leasing period. This offer includes a free LUCIA Equine SAA Kit.

Inquire about leasing options that meet your business requirements:
sales@labmaster.fi.



Labmaster Ltd, based in Kaarina, Finland, is a privately owned company focusing on delivering innovative diagnostics tools for applications in routine clinical diagnostics and research. Labmaster develops and commercializes the next generation of versatile miniaturized detection tools based on its proprietary cathodic electrochemiluminescence-technology (LM-CECL). The superiority of LM-CECL technology allows it to overcome the performance limitations of previous generation's tests relating to sensitivity, accuracy, usability and cost efficiency.

Labmaster's quality management system is ISO 13485:2016 and 9001:2015 certified for the design and manufacture of *in vitro* diagnostic devices for detection of viral infections, immune response and inflammatory status as an aid in diagnosis.

LABMASTER
 ●●● medical devices

Labmaster Ltd
 Rauhalinnantie 31
 20780 Kaarina, Finland
 Tel : +358 22 760 555
 E-mail : sales@labmaster.fi
www.labmaster.fi