

LABMASTER LUCIA™ Veterinary Diagnostics

QUICK AND ACCURATE RESULTS - WITH EASE



LABMASTER LUCIA™ VET ANALYZER

- measures analytes in whole blood
- a quantitative result in a few minutes
- a versatile platform for multiple test kits
- maintenance-free
- easy to transport

POINT-OF-CARE TEST KITS including

- TEST CASSETTE with reagents
- DILUTION TUBE for preparing the sample
- NFC CARD with test-specific parameters

FOAL IgG TEST KIT

- for measuring IgG levels in whole blood to assess failure of transfer of passive immunity (FTPI)

EQUINE SAA TEST KIT

- for measuring the concentration of Equine Serum Amyloid A in whole blood to assess infection and inflammation status in horses

Check out our pipeline veterinary tests at:
www.labmaster.fi

Legal disclaimer

Products designed, developed and manufactured by Labmaster Ltd. meet their specifications when stored and used according to the instructions for use. Labmaster LUCIA™ is a trademark designed and owned by Labmaster Ltd.
Products presented may not be available in all countries. Please check availability from your local representative.

LABMASTER LUCIA™

Foal IgG



INTENDED USE

The Labmaster LUCIA™ Foal IgG test is an *in vitro* veterinary diagnostic test for the quantitative determination of Immunoglobulin G (IgG) from whole blood to assess failure of transfer of passive immunity (FTPI).


KIT CONTENTS (10 tests or 20 tests per Kit)

Kits include test cassettes with reagents, dilution tubes, NFC Card and Instructions for Use. The 10 µL Pipette and 10 µL pipette tip rack are available separately.


SPECIFICATION

SAMPLE TYPE	Whole Blood
SAMPLE VOLUME	10 µL
MEASUREMENT TIME	6 minutes
MEASURING RANGE	100–1500 mg/dL
STORAGE	+2 – +8°C

SIMPLE WORKFLOW OF POINT-OF-CARE TESTING

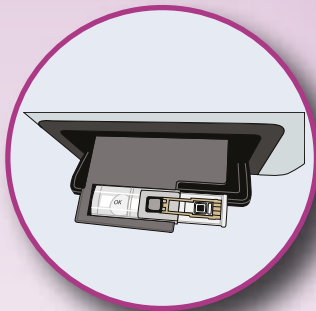

Insert the cassette into the analyzer

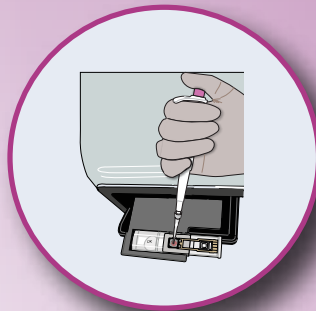

Dispense the prepared sample into the cassette



Result can be transferred to Windows or Android App

Take advantage of the Labmaster LUCIA™ Vet Analyzer leasing programme!

We offer a full ownership of the analyzer after a 24-month leasing period.

Inquire about leasing options that meet your business requirements:
sales@labmaster.fi

Labmaster Ltd, based in Kaarina, Finland, is a privately owned company focusing on delivering innovative diagnostics tools for applications in routine clinical diagnostics and research. Labmaster develops and commercializes the next generation of versatile miniaturized detection tools based on its proprietary cathodic electrochemiluminescence-technology (LM-CECL). The superiority of LM-CECL technology allows it to overcome the performance limitations of previous generation's tests relating to sensitivity, accuracy, usability and cost efficiency.

Labmaster's quality management system is ISO 13485:2016 and 9001:2015 certified for the design and manufacture of *in vitro* diagnostic devices for detection of viral infections, immune response and inflammatory status as an aid in diagnosis.

LABMASTER
●●● medical devices

Labmaster Ltd
Rauhalinnantie 31
20780 Kaarina, Finland
Tel : +358 22 760 555
E-mail : sales@labmaster.fi
www.labmaster.fi